




Curiox C-FREE Pluto MT System

SPECIFICATION SHEET

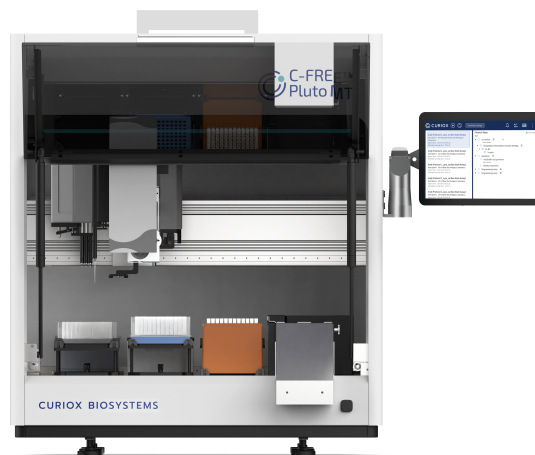

The Curiox C-FREE™ Pluto MT System is an affordable, configurable automation system for cell sample preparation workflows. This automated liquid handling platform enables cell washing and staining processes using standard format 96-well microplates, deep-well plates, tubes and pipette tips. Minimal optimization is needed to port your assay onto the system. The Pluto MT system is a mid-throughput instrument that offers ease of use with Curiox assay protocol scripting services and supports personalized automation, enhancing efficiency, consistency, and reliability.

Performance Specifications

Pipette Precision	5 μ L CV<4%; 200 μ L CV<1%
-------------------	------------------------------------

Maximum Pipetting Volume	1 - 1000 μ L
--------------------------	------------------

Sample Throughput	24 - 72 wells per run
-------------------	-----------------------

Sample Types	Whole blood, bone marrow, mammalian cells, yeast cells
--------------	--

Key Applications	Flow cytometry immunophenotyping assays; cell sorting, multi-omics sample preparation
------------------	--

System Specifications

Dimensions (D x W x H) - System

827 mm* x 790 mm x 895 mm
32 in.* x 31 in. x 35 in.

Dimensions (D x W x H) - Door

690 mm x 470 mm
27.16 in. x 18.5 in.

Weight

140 kg / 309 lbs

Power Supply

100 – 240 VAC, 10A, 50/ 60 Hz

Arm Configurations

1-channel pipetting module
24-channel fixed pipetting module
Gripper module (optional)

Deck Capacity

15 stations and solid waste chute

Pipetting Principle

Air displacement; simultaneous pipetting of 24 channels

Liquid Level Detection

Pressure liquid level detection

Calibration

Engineer calibrated during installation and annual preventive maintenance**

Touchscreen

12 in.

Communication Interface

USB, ethernet

System Configuration

Stand-alone

* Dimensions include waste chute

** Annual Preventive Maintenance is highly recommended to ensure that the system is functioning at peak performance.

Description

Part No.

24 C-FREE, Pluto MT System, 24-channel Liquid Handling Head

CF-PLU-LH08MT-01

See the future of cell sample preparation automation at www.curiox.com/pluto



For Research Use Only. Not for use in diagnostic procedures.

© 2024 Curiox Biosystems, Inc. All rights reserved. C-FREE Pluto is a trademark of Curiox Biosystems, Inc.
09/2024